


SYSTEM ADR III-7 Reciprocating Cylinder/ Reciprocating Holder Method

Simulate the passage of your dose through the GIT by exposing the preparation to changes in media such as pH automatically and through extended longitudinal studies.

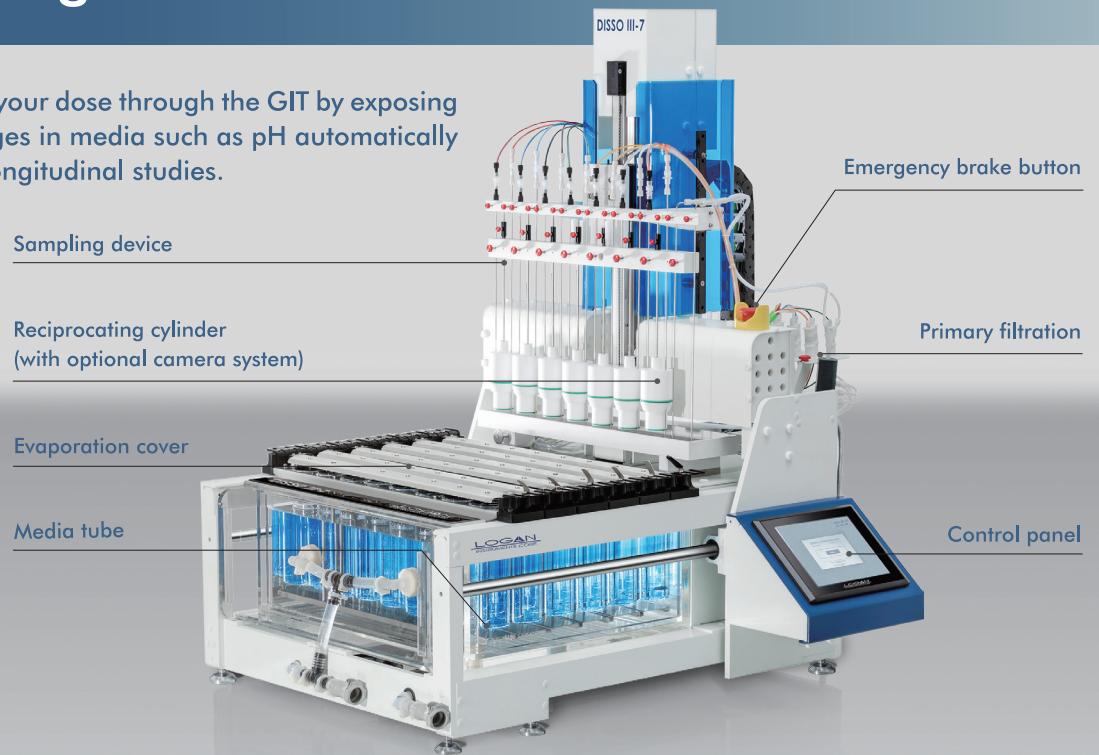

Reciprocating Cylinder/Holder

Reciprocating cylinders/holders of various specifications are optional.



Samples Taken and Media Replacement

Individual sampling probe and individual media replacement probe work precisely to take samples.



Evaporation Cover

Each row of the media tubes is covered.



Camera System

An optional camera system can be provided to observe and record the whole process of drug release.

Auto Cleaning

The sampling and media replacement tubing can be automatically cleaned with the program.

Easy Operation

The standard reciprocating cylinder and all sampling holders are easy to install.



Sampling

Pre-rinse sample lines before taking samples.



Automatic Sampling

DSC-800 dissolution system controller, SCR - DL auto sampler perform automatic sampling.



Customized service

LOGAN not only has USP Apparatus 1-7 dissolution testers, but also provides various customized design services to meet the research & development needs of various new drugs, such as testing of contraceptive rings, cardiac stents, hip stems, and novel nanoparticle formulations.

